



A Dräger and Siemens Company

# Field Service Procedure

Part Number: SP00286



Rev: -

Date: 15 September 2003


© 2003 Draeger Medical, Inc.

## Double Vapor Parking Installation

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)  
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

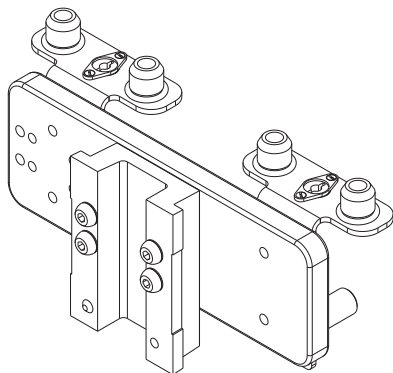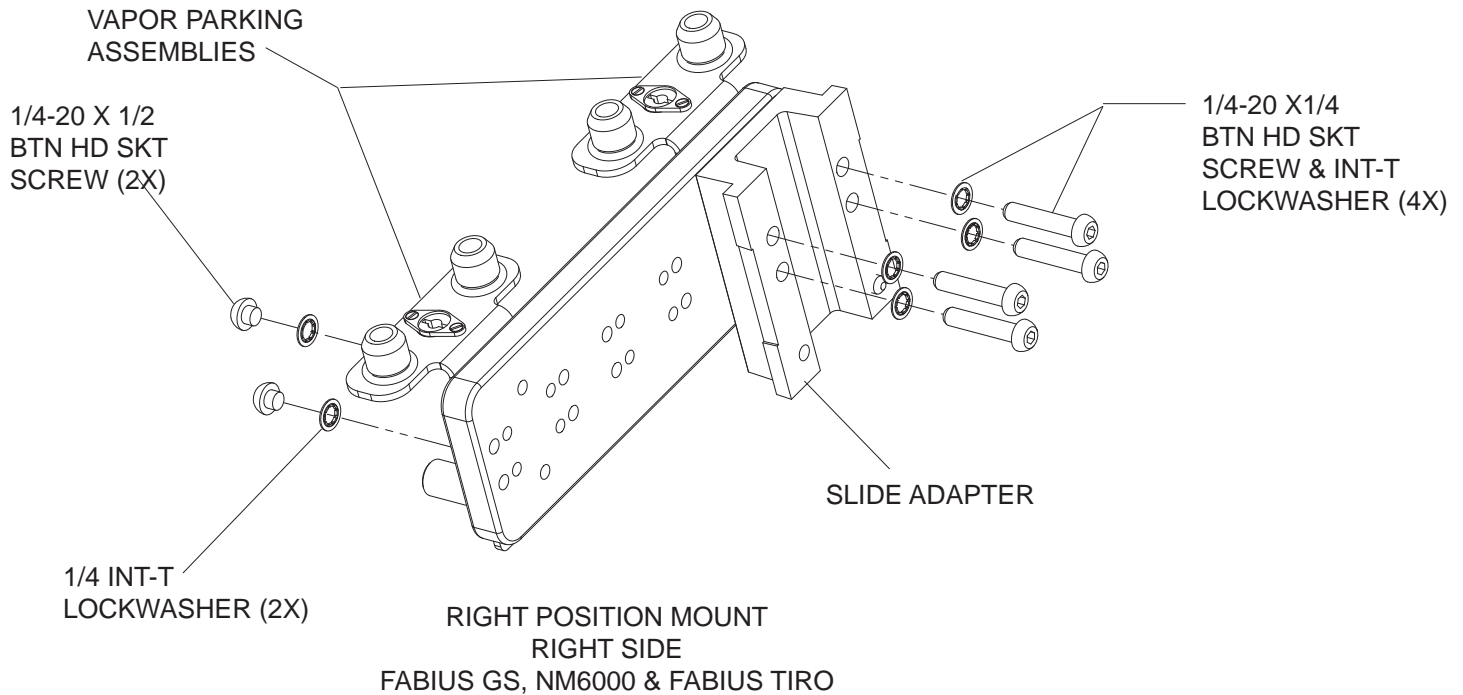

CENTER POSITION MOUNT  
FABIUS GS & NM6000

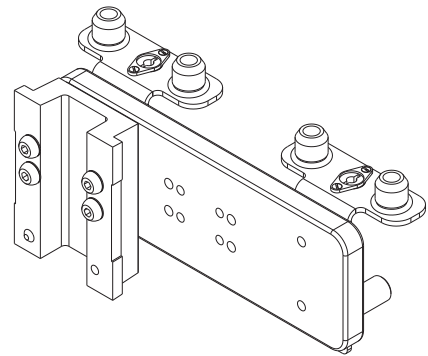

LEFT POSITION MOUNT  
RIGHT SIDE  
FABIUS GS, NM6000 & FABIUS TIRO

1. Select a permitted double vapor parking arrangement (left, center or right position).
2. Attach the slide adapter to the mounting plate with four 1/4-20 x 1-1/4 in. btn hd skt screws and 1/4 int-t lock washers as shown in the illustration.
3. Attach each vapor parking assembly to the mounting plate as shown in the illustration with two 1/4-20 x 1/2 in. btn hd skt screws and 1/4 int-lock washers.


# Installation:Continued:

Note: Machine must be equipped with a GCX rail.


4. Install the double vapor parking accessory in the GCX rail. Position the stop piece in the GCX rail at the selected position and tighten the setscrews.

**WARNING:** Do not install the double vapor parking accessory on the left side GCX rail on a Narkomed 6000 Series Machine.

## TYPICAL NM6000 INSTALLATIONS:



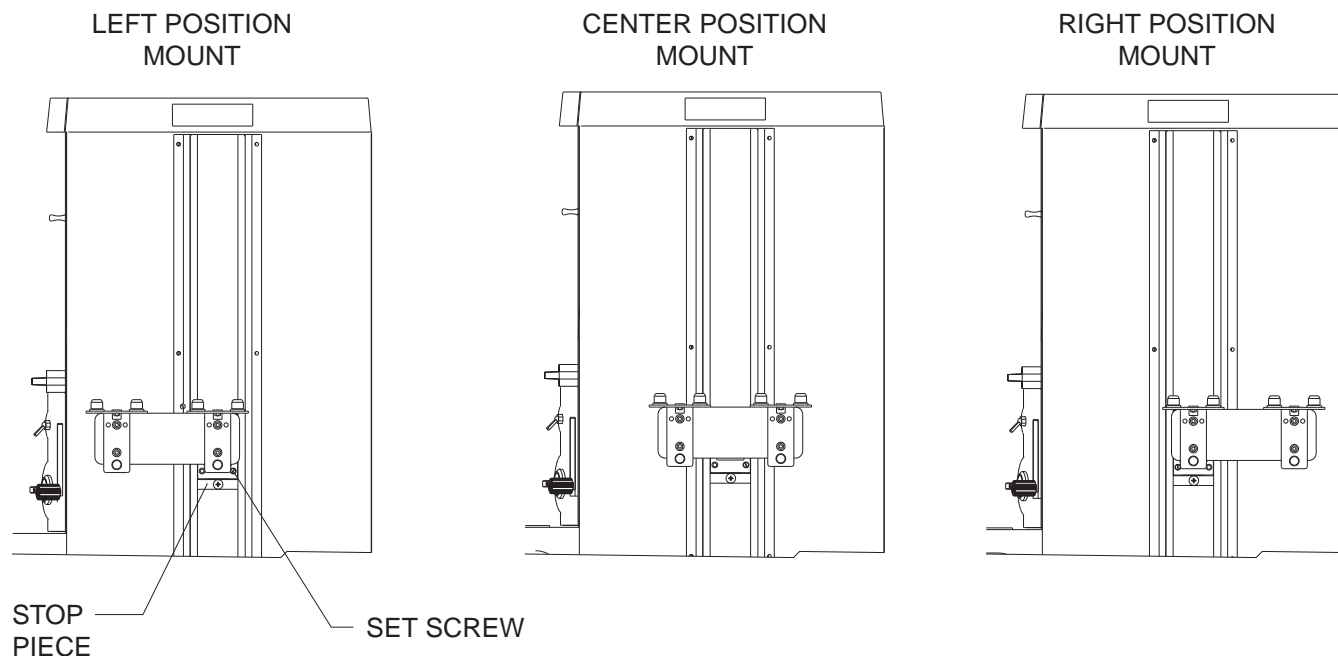


## TYPICAL FABIUS TIRO INSTALLATION:

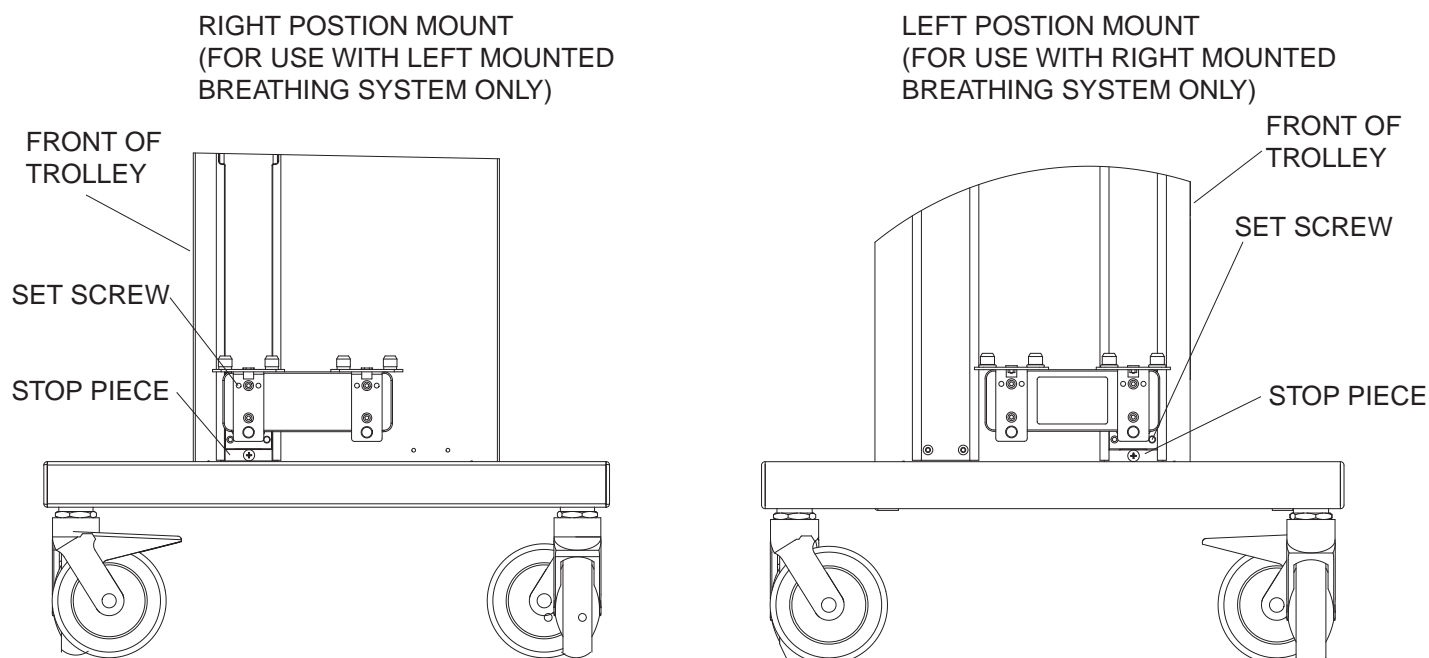

VAPOR PARKING ASM

TYPICAL FABIUS GS INSTALLATIONS:

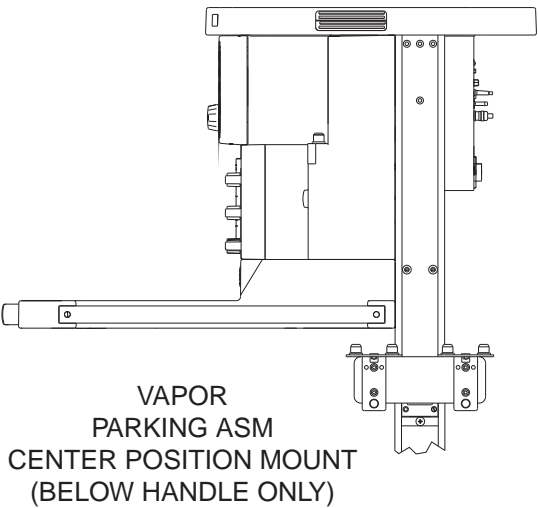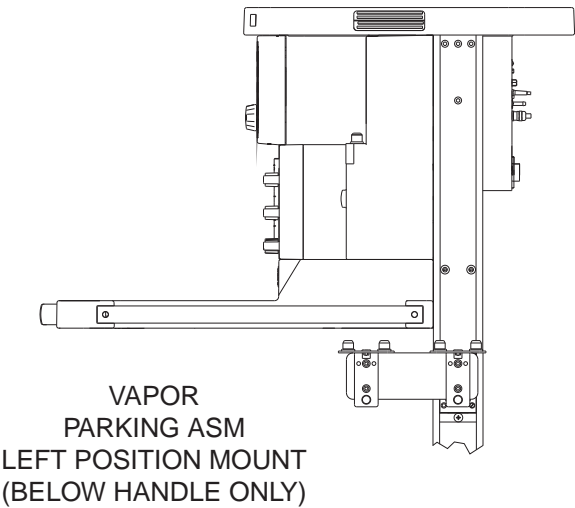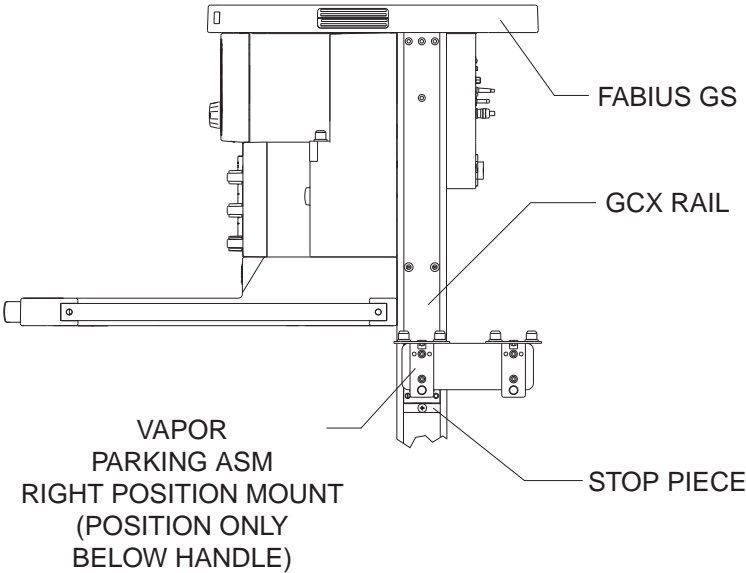

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)  
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)  
[RETURN TO CD-ROM TABLE OF CONTENTS](#)



DrägerService is a division of  
Draeger Medical, Inc.  
3122 Commerce Drive  
Telford, PA 18969  
Tel: (215) 721-5402  
(800) 543-5047  
Fax: (215) 721-5784  
Web: [www.draegermedical.com](http://www.draegermedical.com)  
Printed in the U.S.A.